

Certificate of Analysis

Product Code:	LIS-0480S
Product:	Solus <i>Listeria</i> ELISA (5 plate)
LOT/Batch:	EL094
Date of Manufacture/Test:	18.12.2017
Expiry Date:	17.12.2018

This batch has been quality control tested to the following criteria:

Batch homogeneity: CV%	< 10%.
A panel of 9 non- <i>Listeria</i> sp.	<u>Satisfactory</u>
A panel of 10 <i>Listeria</i> sp.	<u>Satisfactory</u>

Performance of In house Reference Standards

Sample (cfu/ml)	Net Signal OD	OD Specification
D9 <i>L. monocytogenes</i> LMG 16780 10 ⁶	1.420	>0.200
D9 <i>L. monocytogenes</i> LMG 16780 10 ⁷	>2.500	>1.500
C10 <i>L. monocytogenes</i> ATCC 13932 10 ⁶	0.259	>0.200
C10 <i>L. monocytogenes</i> ATCC 13932 10 ⁷	1.693	>1.500
D8 <i>L. monocytogenes</i> ATCC 19114 10 ⁶	1.037	>0.200
D8 <i>L. monocytogenes</i> ATCC 19114 10 ⁷	>2.500	>1.500
D7 <i>L. monocytogenes</i> ATCC 19111 10 ⁶	0.210	>0.200
D7 <i>L. monocytogenes</i> ATCC 19111 10 ⁷	1.693	>1.500
G3 <i>L. innocua</i> (wild type) 10 ⁶	0.635	>0.200
G3 <i>L. innocua</i> (wild type) 10 ⁷	>2.500	>1.500

Solus Scientific Solutions Ltd

9 Mansfield Networkcentre, Millennium Business Park, Concorde Way, Mansfield, Nottinghamshire NG19 7JZ

Tel +44 (0)1623 429701 Fax +44 (0)1623 620977 Email enquiries@solusscientific.com

Web www.solusscientific.com


Registered Address: East Kilbride Laboratory, Torus Building, Rankine Avenue, East Kilbride, Glasgow, Lanarkshire G75 0QF Company Reg. No: SC361616, Registered in Scotland VAT No 982 2304 24

Net Optical Densities of samples grown in RELM at 30.0°C for 24 hours and boiled at 85-100°C for 20 minutes.

	1	2	3
A	Negative control 0.049	<i>D5 L. ivanovii</i> >2.500	<i>E9 Bacillus cereus</i> 0.046
B	Positive control >2.500	<i>D9 L. monocytogenes</i> >2.500	<i>A8 E. coli</i> 0.037
C	Substrate blank 0.061	<i>D8 L. monocytogenes</i> >2.500	<i>C1 Staphylococcus aureus</i> 0.064
D	<i>C10 L. monocytogenes</i> >2.500	<i>D7 L. monocytogenes</i> >2.500	<i>E2 E. coli (WT)</i> 0.041
E	<i>G10 L. innocua</i> >2.500	<i>D10 L. seeligeri</i> >2.500	<i>E8 Brochothrix thermosphacta</i> 0.045
F	<i>D2 L. grayii (BHIB)</i> >2.500	<i>E7 Kurthia zopfii</i> 0.050	<i>E6 Enterococcus faecalis</i> 0.056
G	<i>D1 L. monocytogenes</i> >2.500	<i>F1 Streptococcus bovis</i> 0.043	Absorbances read at 450nm
H	<i>D6 L. welshimeri</i> >2.500	<i>F10 Enterococcus faecalis</i> (WT) 0.061	

Production Batch is approved for release as it gives satisfactory QC performance.

Signed: ~~Operations director/Technical director/CEO~~

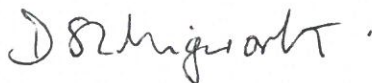

Date: 22/12/17.

Solus Scientific Solutions Ltd

9 Mansfield Networkcentre, Millennium Business Park, Concorde Way, Mansfield, Nottinghamshire NG19 7JZ
 Tel +44 (0)1623 429701 Fax +44 (0)1623 620977 Email enquiries@solusscientific.com
 Web www.solusscientific.com

Registered Address: East Kilbride Laboratory, Torus Building, Rankine Avenue, East Kilbride, Glasgow,
 Lanarkshire G75 0QF Company Reg. No. SC361616, Registered in Scotland VAT No 982 2304 24