

Certificate of Analysis

Product Code:	LIS-0480S
Product:	Solus <i>Listeria</i> ELISA (5 plate)
LOT/Batch:	EL090
Date of Manufacture/Test:	16.10.2017
Expiry Date:	15.10.2018

This batch has been quality control tested to the following criteria:

Batch homogeneity: CV%	< 10%.
A panel of 9 non- <i>Listeria</i> sp.	<u>Satisfactory</u>
A panel of 10 <i>Listeria</i> sp.	<u>Satisfactory</u>

Performance of In house Reference Standards

Sample (cfu/ml)	Net Signal OD	OD Specification
D9 <i>L. monocytogenes</i> LMG 16780 10 ⁶	0.585	>0.200
D9 <i>L. monocytogenes</i> LMG 16780 10 ⁷	>2.500	>1.500
C10 <i>L. monocytogenes</i> ATCC 13932 10 ⁶	0.237	>0.200
C10 <i>L. monocytogenes</i> ATCC 13932 10 ⁷	1.538	>1.500
D8 <i>L. monocytogenes</i> ATCC 19114 10 ⁶	0.727	>0.200
D8 <i>L. monocytogenes</i> ATCC 19114 10 ⁷	>2.500	>1.500
D7 <i>L. monocytogenes</i> ATCC 19111 10 ⁶	0.221	>0.200
D7 <i>L. monocytogenes</i> ATCC 19111 10 ⁷	1.629	>1.500
G3 <i>L. innocua</i> (wild type) 10 ⁶	0.278	>0.200
G3 <i>L. innocua</i> (wild type) 10 ⁷	1.691	>1.500

Solus Scientific Solutions Ltd

9 Mansfield Networkcentre, Millennium Business Park, Concorde Way, Mansfield, Nottinghamshire NG19 7JZ

Tel +44 (0)1623 429701 **Fax** +44 (0)1623 620977 **Email** enquiries@solusscientific.com

Web www.solusscientific.com


Registered Address: East Kilbride Laboratory, Torus Building, Rankine Avenue, East Kilbride, Glasgow, Lanarkshire G75 0QF Company Reg. No. SC361616, Registered in Scotland VAT No 982 2304 24

Net Optical Densities of samples grown in RELM at 30.0°C for 24 hours and boiled at 85-100°C for 20 minutes.

	1	2	3
A	Negative control 0.049	<i>D5 L. ivanovii</i> >2.500	<i>E9 Bacillus cereus</i> 0.037
B	Positive control >2.500	<i>D9 L. monocytogenes</i> >2.500	<i>A8 E. coli</i> 0.035
C	Substrate blank 0.052	<i>D8 L. monocytogenes</i> >2.500	<i>C1 Staphylococcus aureus</i> 0.041
D	<i>C10 L. monocytogenes</i> >2.500	<i>D7 L. monocytogenes</i> >2.500	<i>E2 E. coli (WT)</i> 0.036
E	<i>G10 L. innocua</i> >2.500	<i>D10 L. seeligeri</i> >2.500	<i>E8 Brochothrix thermosphacta</i> 0.057
F	<i>D2 L. grayii (BHIB)</i> 2.204	<i>E7 Kurthia zopfii</i> 0.042	<i>E6 Enterococcus faecalis</i> 0.053
G	<i>D1 L. monocytogenes</i> >2.500	<i>F1 Streptococcus bovis</i> 0.037	Absorbances read at 450nm
H	<i>D6 L. welshimeri</i> >2.500	<i>F10 Enterococcus faecalis (WT)</i> 0.054	

Production Batch is approved for release as it gives satisfactory QC performance.

Signed: Operations director/Technical director/CEO

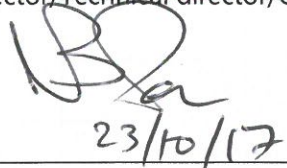

Date:

23/10/17

Solus Scientific Solutions Ltd

9 Mansfield Networkcentre, Millennium Business Park, Concorde Way, Mansfield, Nottinghamshire NG19 7JZ
Tel +44 (0)1623 429701 **Fax** +44 (0)1623 620977 **Email** enquiries@solusscientific.com
Web www.solusscientific.com

Registered Address: East Kilbride Laboratory, Torus Building, Rankine Avenue, East Kilbride, Glasgow, Lanarkshire G75 0QF Company Reg. No. SC361616, Registered in Scotland VAT No 982 2304 24